


# Instructions for the Pre/Post Event Monitor

## ATTACHING THE WIRES:

1. Place electrodes on clean, dry skin as shown in diagram.
2. Attach **WHITE** wire to right side of chest below collarbone
3. Attach **RED** wire under left breast on ribcage

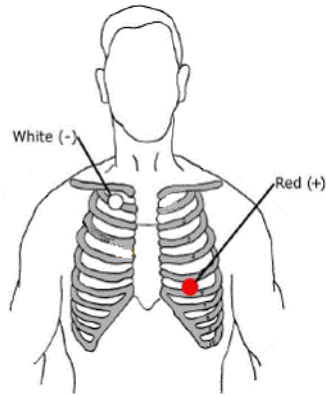

## MONITOR SET- UP:

1. Connect the round end of the patient cable to the monitor.
2. Three quick beeps will be heard indicating a good connection.

## WHEN SYMPTOMS OCCUR:

1. Press **RECORD** button for approximately 2 seconds until a short beep is heard.
  - If no beep is heard, mute is turned on.
2. Relax and sit or stand as still as possible while recording for 30 seconds.
3. When recording stops, you will hear a telephone ring sound.

## TRANSMITTING:

1. Set the Monitor on a flat surface with the record button facing up.
2. Call 1-800-747-4455 from a land line phone. A technician will help you through the process.
3. When instructed, push the **Send** button for two seconds.
4. Pressing the **Send** button for two seconds during a transmission will abort the transmission.



## TO MUTE SOUND:

Mute will silence the recorder except when recording events

1. Press the mute button for two seconds to change the mute setting
  - A single beep will sound when mute is turned on.
  - A double beep will sound when the mute is turned off.
  - Mute default is off when a battery is inserted
  - Mute defaults off after sending events

## RESETTING THE MONITOR:

Push the **SEND** and **RECORD** buttons simultaneously for three seconds.

- The monitor will beep once when all events are erased.
- The monitor will also automatically erase events when the **RECORD** button is pushed for two seconds after all events have been sent successfully