


Instructions for the 180-60 Multi-Event Monitor

ATTACHING THE WIRES:

Place Electrodes on clean, dry skin in location shown in diagram. Attach WHITE wire(-) to right side of chest.

Attach GREEN wire (ground) to left side of chest.

Attach BLACK wire (+) under left breast on ribcage.

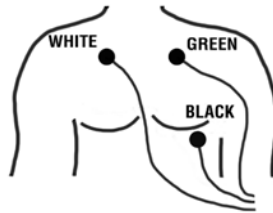

Follow Color Scheme not markings on wires



MONITOR SET-UP:

1. Connect lead wires to monitor by inserting the connector into the single port on top of the monitor.
2. Install battery. Unit will beep 5 times.
3. To erase memory, see **RESETTING THE MONITOR** section. This will clear the monitor and ready it for recording an event.

WHEN SYMPTONS OCCUR

1. Press BLUE Record button until a beeping sound is heard.
2. Relax and sit or stand as still as possible while recording for 60 seconds.
3. When recording stops, you will hear the beeping sound again and the GREEN light will flash for approximately 30 seconds.
4. Only one event can be recorded on this device.

TRANSMITTING:

1. Press the GREEN send button on side of monitor until you hear a constant tone, then release button and hold phone in position until tone stops.
2. After tone stops, monitor will continue beeping and green light will flash until monitor is reset.

RESETTING THE MONITOR:

1. RESET monitor by holding both BLUE and GREEN buttons down at the same time.
2. HOLD until you hear 5 beeps. Monitor is now reset and ready to record again.