

Instructions for the Multi-Event Monitor


ATTACHING THE WIRES:

Place Electrodes on clean, dry skin in location shown in diagram.

Attach WHITE wire (-) to right side of chest.

Attach GREEN wire (ground) to left side of chest.

Attach BLACK wire (+) under left breast on ribcage.

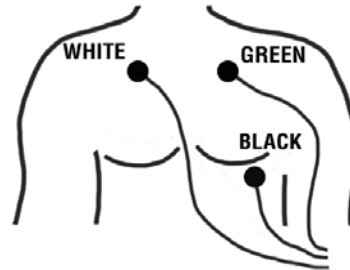

MONITOR SET-UP:

Connect lead wires to monitor by inserting the connector into the single port on top of the monitor.

Install battery. Unit will beep 5 times.

To erase memory, see **RESETTING THE MONITOR** section. This will clear the monitor and ready it for recording an event.

WHEN SYMPTONS OCCUR

Press BLUE Record button until a beeping sound is heard.

Relax and sit or stand as still as possible while recording for 30 seconds.

When recording stops, you will hear the beeping sound again and the GREEN light will flash for approximately 30 seconds.

Wait 45 seconds before recording another event.

Number of flashes equals number of events stored – One Flash = 1 stored event; Two flashes – 2 events; Three flashes = 3 events; Four flashes = 4 events

TRANSMITTING:

Place Monitor on a flat Surface.

Place phone mouthpiece over hole on front of monitor as instructed by technician.

Press the GREEN send button on side of monitor until you hear a constant tone, then release button and hold phone in position until tone stops.

Transmission takes approximately one minute per event.

After tone stops, pick up the phone and listen to technician for further instructions. Monitor will continue beeping and green light will flash until monitor is reset.

RESETTING THE MONITOR:

RESET monitor by holding both BLUE and GREEN buttons down at the same time.

HOLD until you hear 5 beeps. Monitor is now reset and ready to record again.

Wait 45 seconds before recording the next event.