


Instructions for the Genesis Multi-Memory Event

ATTACHING THE WIRES:

1. Place Electrodes on clean, dry skin as shown in diagram.
2. Your physician may prefer the 2 or 3 wire configuration.
3. For the 2 wire use the red and white position only

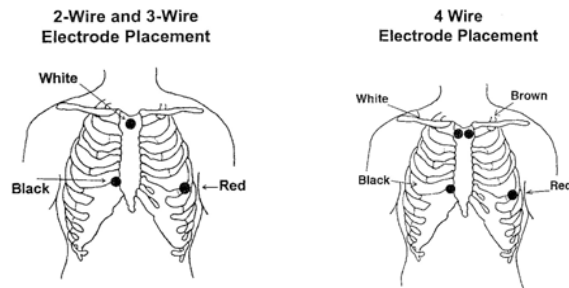

MONITOR SET-UP:

1. Connect lead wires to monitor by inserting the connector into the single port on top of the monitor.
2. Install battery. Unit will beep 5 times.
3. To erase memory, see **RESETTING THE MONITOR** section. This will clear the monitor and ready it for recording an event.

WHEN SYMPTONS OCCUR

1. Press **RED STORE** button until a conformation beep is heard.
2. Relax and sit or stand as still as possible while recording.
3. When recording stops, the LCD will display how many events have been stored.
4. Wait 45 seconds before recording another event.



TRANSMITTING:

1. Press the **SEND** button, then release button and hold phone in position until tone stops.
2. Transmission takes approximately one minute per event, or can be increased by doctor's request.
3. During transmission monitor will display "Sending" on the display. The transmission is complete when the tones stop and the display reads "Done Sending."

RESETTING THE MONITOR:

1. Events must be transmitted in order to reset the monitor.
2. **RESET** monitor by holding the **STORE** button down for 3 seconds.
3. **HOLD** until you hear a series of beeps. Monitor is now reset and ready to record again.
4. Wait 45 seconds before recording the next event.